

TOP GLOVE SDN. BHD

WORK INSTRUCTION
Finished Product Testing : Aging Test

Purposes/Function/Objective

- To observe the glove quality either stickiness or discolor for powder free and nitrile glove.

Materials/ Chemicals/

Tools/Equipment

1. Inner (sample gloves)
2. Oven 70 ° C

Form/s:

- LA/F16

References :

- nil


Procedures



Take few pieces of gloves or 1 inner (100 pieces) sample from any line or customer brand.



Put inside glove wallet or inner box

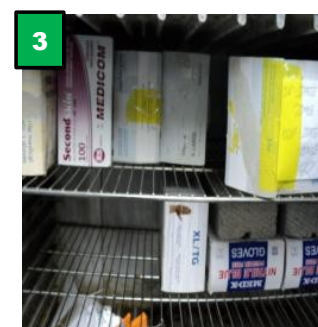

Keep inside oven for 22 hours at 70 °C

5
Let the sample cool to room temperature and do the observation.

Prepared by:	Verified by:	Approved by:
Noor Azura Binti Azman QA/MQC Asst Supervisor	Al-Fadilah Mohamad QA Executive II	Noor Akilah Saidin QA Deputy General Manager

(Effective Date : 02/05/2011)

