TOP GLOVE SDN. BHD

WORK INSTRUCTION

Finished Product Testing: Tensile test

Purposes/Function/Objective

Procedures

- To check physical properties of glove

<u>Materials/ Chemicals/</u> <u>Tools/Equipment</u>

- 1. dumb bell cutter
- 2. Cutter machine
- 3. Thickness gauge
- 4. Tensile machine

Specification

- 1. There are many types of dumb bell cutter(ISO, ASTM, EN, ABNT)
- 2.Tensile machine test will automatically reset to its original position when the specimen breaks.
- 3. For the next sample, repeat procedure no 4,5 & 7
- Click "report" at the top panel when finish testing.
 Adjust title and other relevant info inside the report.
- 5. Save and add file name. Print report if needed.

Form/s:

- NIL

References:

- nil

(Effective Date: 02/05/2011)




Cut glove with dumb bell cutter at palm area follow by standard.



Measure thickness of the specimen at 3 position within 2.5cm marks.



Put glove inside the cutter machine and press.

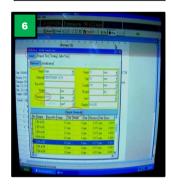

Open the tensile machine file at the computer.



Separate each piece of dumb bell specimen and label each specimen.



Draw 2 bench marks of 2.5cm apart of center narrow portion of the specimen.







Fix the dumb bell sample at tensile machine. The sample are moving upward pulling along the specimen.

Prepared by:	Verified by:	Approved by:
Noor Azura Binti Azman	Al-Fadilah Mohamad	Noor Akilah Saidin
QA/MQC Asst Supervisor	QA Executive II	QA Deputy General Manager