

TOP GLOVE SDN. BHD

WORK INSTRUCTION

Finished Product Testing : Protein Test Method

Sample Preparation (ASTM D 5712)

Purposes/Function/Objective

To determine protein content in

sample glove

Materials/ Chemicals/

Tools/Equipment

-Scissors

- glove sample (3 pieces)

Reference


ASTM D 5712

Form/s

-LA/F22

(Effective Date : 02/05/2011)

Cut gloves into small pieces (finger-10, palm/cuff-24) as per below picture:

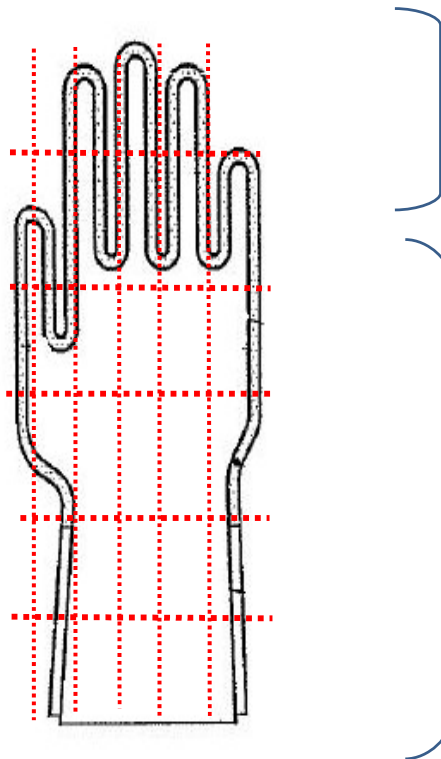

Part 1 : 10 pcs

Part 2 : 12 x 2 (2sides)  
= 24 pcs

**TOTAL =34 pieces**

Prepared by:	Verified by:	Approved by:
Noor Azura Binti Azman QA/MQC Asst Supervisor	Al-Fadilah Mohamad QA Executive II	Noor Akilah Saidin QA Deputy General Manager

